

SMALL ANIMAL EQUIPMENT

Y1394 BONART MICRO-MOTOR / POLISHER

The Micro-motor/polisher makes up an essential part of the scaling program. Polishing (with prophy paste and prophy cups) smooths the tooth surface and helps prevent plaque building up again too rapidly after scaling (scaling leaves the tooth surface rough). The micro-motor is a low speed, high torque motor, with variable speeds from 300 to 30,000 RPM. Lower speeds up to 2,500 RPM are used for polishing, while higher speeds are used for cutting.

Y1394 includes the following:

- Y1380 NSK AR-Y latch-type head for burs
- Y1391 NSK EC handpiece for micro-motor (no heads)
- Y1399 NSK NAC-Y snap-on prophy head for polishing
- Y1382 Prophylactic cups, rubber (10's)
- Y1383 Prophylactic paste (10's)
- Y1384 Tungsten carbide round burr #3
- Y1384/701 Cutting Burr #701 (small)
- Y1384/703 Cutting Burr #703 (large)
- Y1391/HP NSK HP Straight hand piece (optional)
- Y1384/HP023 HP Surgical round burr (optional)
- Y1384/HP016 HP Surgical cutting burr (optional)
- Y1375 Lubricant oil for NSK handpieces (480ml)
- Y1376 E-Type nozzle to use with lubricant oil



Y1394



Y1384 (round) • Y1384/701 & 703 (s/cutting)



Y1380 • Y1391 • Y1393



Y1382 • Y1383



Y1391/HP • Y1384/HP023 • Y1384/HP016



Y1375



Y1376

Y1390 BONART ULTRASONIC SCALER

This ultrasonic scaler incorporates the latest in ultrasonic technology at an affordable price. Has adjustable water and power control that provides faster cleaning in any situation.

Y1390 includes the following:

Y1395 Narrow/universal tip scaling insert

Y1396 Spatula/beavertail tip scaling insert

Y1392 Female water connector

Y1392/m Male water connector
water pump/bottle

Y1393 BONART ULTRASONIC SCALER & MICROMOTOR / POLISHER COMBINATION

Combines all the features and items of the scaler (Y1390) and micro-motor/polisher (Y1394) in one durable compact unit to give versatility and scope for successful veterinary dentistry.

See page 137 for cleaning, oiling and general maintenance of NSK handpieces and heads.

Manufacturing Warranty Against Defects


- *Equipment:* 1 year from date of purchase and excludes expendable goods (e.g. inserts, cable, foot switch, fuses, etc).
- *NSK handpieces and heads:* 6 months from date of purchase.

Please note that the warranty will be forfeited if the units are opened by anyone other than Kyron's authorised technician in Johannesburg.

Water pressure to the units needs to be regulated — maximum water pressure should not exceed 30psi.



Y1390

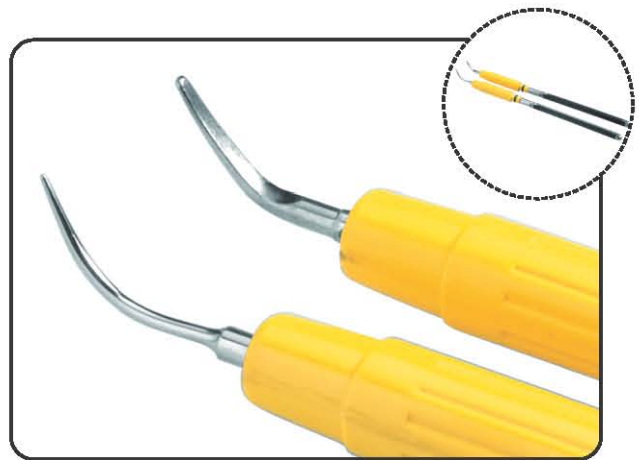

Y1395 • Y1396



Y1393

FACE SHIELD

Head band with plastic visor to protect the operator from inhaling the fine H₂O mist (filled with potentially harmful bacteria) created when scaling.



Face Shield

DENTAL EXTRACTOR


Y1316 CURVED: For lower molars
Length: 16,5cm

Y1317 FRAGMENT: For tooth fragments and small teeth.
Length: 16,5cm

Y1315 STRAIGHT: For upper molars.
Length: 16,5cm



Y1316

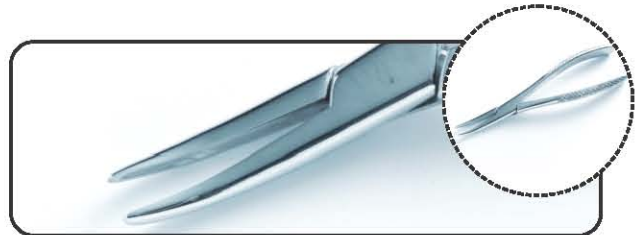

Y1317

DENTAL SCALER


Y1322 Sickle Shape
Stainless steel – with one sharp sickle end.
Length: 15,5cm

Y1321 Universal Double Ended
Stainless steel – with two flat-edged curved ends.
Length: 17cm

Y1320 Universal Single Ended
Stainless steel – with one flat-edged curved end.
Length: 15,5cm



Y1315

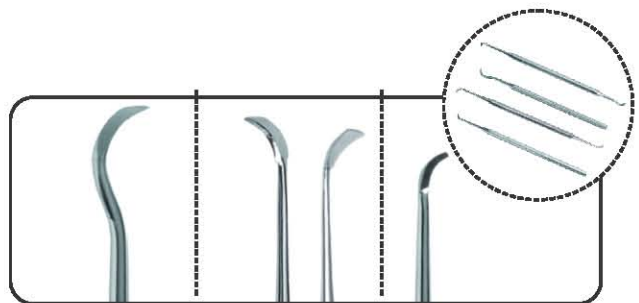

Y1322 • Y1321 • Y1320

Y1323 Calculus Remover
For the removal of very hard deposits of calculus which often develops on the carnassial and other teeth.



Y1323

SECTION D – DENTAL


DENTAL ROOT ELEVATOR

Length: 14cm

Size:

Y1309 2mm

Y1310 3mm

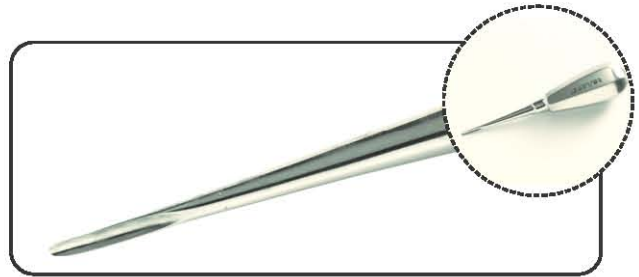

Y1309


WINGED DENTAL ELEVATOR – (Luxator)

Length: 14cm

Size:

Y1311 2mm

Y1312 4mm

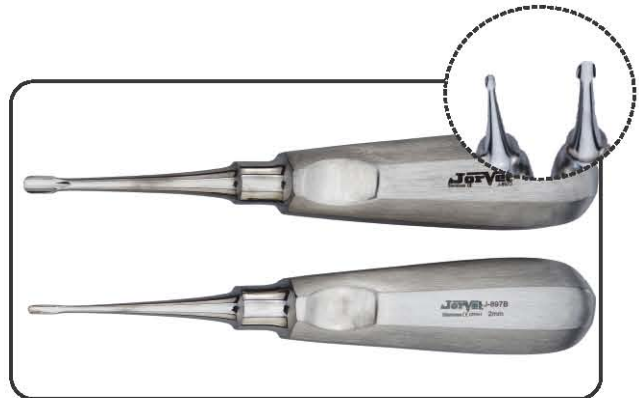

Y1311 • Y1312

ORAL HYGIENE

Y1300 ORAL HYGIENE KIT


Direct brushing of a dog's teeth and gums helps prevent such health problems as doggie breath, loose or infected teeth and gingivitis. Contains a finger-brush and a gum massager.



Y1300

Y1302 PET DENT® FINGER BRUSH


Convenient finger-fitting toothbrush in a single pack. Use with PET DENT® Toothpaste (sold separately).

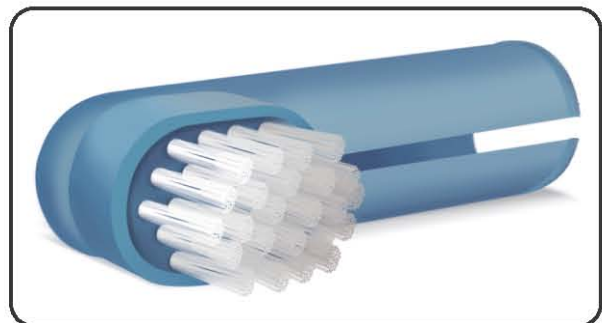

Y1302

Y1301 PET DENT® TOOTHBRUSH

This unique “double ended” toothbrush with its special long handle (22cm) will reach far back in the dog's mouth. Use with PET DENT® Toothpaste (sold separately).

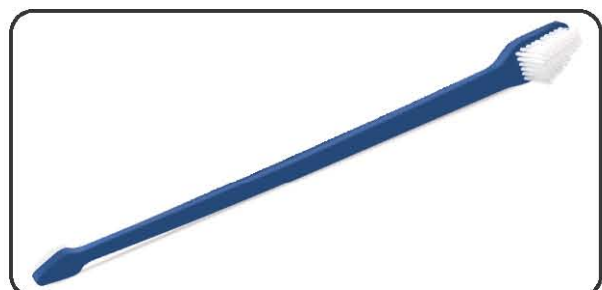

Y1301