


### CLINICAL

- Y2505 PATELLA HAMMER, TAYLOR'S  
A stainless steel reflex hammer with rubber head.  
Size: 20cm

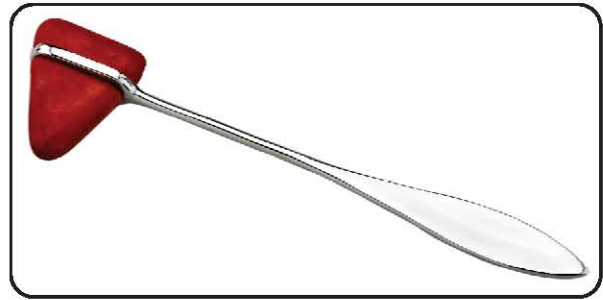

Y2505

### STETHOSCOPES

#### OESOPHAGEAL STETHOSCOPE

A flexible plastic tube with a hollow core used to monitor respiratory and heart sounds when access to the chest is difficult (flexibility of the plastic means traumatic damage is unlikely).

The distal end inserts into the oesophagus and contains a number of small holes with a thin plastic covering (to prevent fluids entering the tube), while the proximal end (fitted with a luer mount) transmits the respiratory and heart sounds. The luer mount fits onto the end of a normal stethoscope (*not included*) (simply remove the stethoscope head and insert the luer mount into the tube). Tube length: 50cm

Sizes:

- Y2561 12FG (2,7mm id x 4mm od)  
Y2563 18FG (4,1mm id x 6mm od)  
Y2564 24FG (5,5mm id x 8mm od)



Y2561 • Y2563 • Y2564

**Y2599 SPRAGUE-RAPPAPORT  
STETHOSCOPE**

This extremely versatile stethoscope is joint-less, dual-tubed and dual-sided — the large side of the chestpiece is 45mm in diameter while the small side of the chestpiece is 30mm in diameter (switch easily between the sides by rotating the head 180°).

Tube length: 48cm

*The stethoscope comes with two extra diaphragms (large and small), three bell-type chestpieces as well as three sets of spare ear tips in different sizes.*



**Y2599**

**Y2562 ULTRASCOPE STETHOSCOPE**

The UltraScope is a pressure-sensitive single-sided single-tube stethoscope, which offers efficiency in auscultation with superior acoustics (use very light pressure for low frequency and firmer pressure for higher frequency responses). The unique, rounded, acrylic head (45mm in diameter) is also crack resistant. Latex-free, cardiology style tubing available in various colours and patterns (always available in paw-print & tiger eye designs, but other designs are available on import – see [www.ultrascope.com](http://www.ultrascope.com). This option takes a few weeks).

*The stethoscope comes with an extra set of ear tips and a lifetime warranty for the original owner*



**Y2562**

**Y2560 STANDARD DOCTOR'S STETHOSCOPE**

Joint-less, single-tube stethoscope – the diaphragm side of the chestpiece is 33mm in diameter while the bell side of the chestpiece is 25mm in diameter. The bell comes fitted with a chill ring.

*This stethoscope comes with an extra set of ear pieces and a spare diaphragm*



**Y2560**

**THERMOMETERS**

**ENVIRONMENTAL**

**Y2570** Alcohol thermometer (-10°C - 110°C)

**RECTAL**


**Y2572** Digital, waterproof thermometer

**Y2571** Digi-Temp thermometer with fast accurate readings on LED display

**Y2579** KYRON DIGITAL THERMOMETER

This professional veterinary thermometer provides instant temperature measurement for all animals — simple, accurate, dependable measurements every time. This fast-testing digital thermometer provides a temperature reading in 8-20 seconds with an accuracy of approximately 0.1°C (in the range 35°C - 39°C) to 0.2°C (<35°C or >39°C). An ergonomically-designed instrument, the dimensions are 190mm H x 110mm probe length x 35mm W.

**Y2573** Standard veterinary thermometer (mercury 37°C - 43°C)

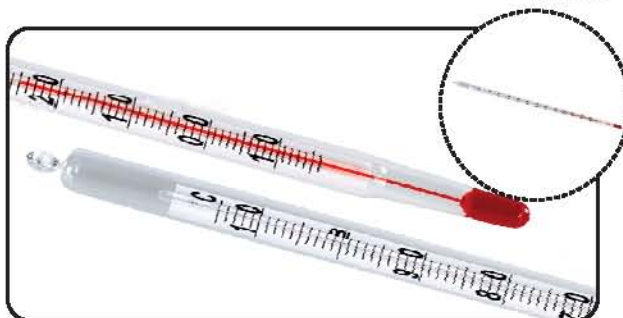

**Y2570**



**Y2572**



**Y2571**



**Y2579**



**Y2573**



**Y6300 TONOMETER, SHIOTZ**  
 Provides accurate measurement of intra-ocular pressure to help with the diagnosis of e.g. glaucoma, uveitis. For your convenience, a table for measuring pressure has been included in this catalogue — see page 130.



**Y6300**

**Y2500 WOODS LAMP**  
 Ultraviolet (black light) light used in the diagnosis of certain fungal or ringworm (*Microsporum canis*) infections that will fluoresce green. Small and portable; uses AA batteries.



**Y2500**

**Y2501 WOODS LAMP REPLACEMENT BULB**

## LABORATORY GENERAL

**Y2275 CLINICAL REFRACTOMETER**  
 A high-quality clinical hand-held refractometer for quick and accurate urine specific gravity and serum protein measurements. Adjustable focus and comes in a protective case.



**Y2275**

**Y2270 COPLIN JARS**  
 Polyethylene jar for holding 6-8 standard slides at a time during staining procedures.

*Also available: Kyron Stains Range which includes Kyroquick set, Sternheimer Malbin, Methylene Blue, etc.*


## FAECAL LOOPS

Inexpensive plastic loop with slotted ends for quick hygienic faecal samples. Sizes:

**Y2530** 24cm (small, medium & large dogs)  
**Y2531** 11cm (very small dogs & puppies)



**Y2270**

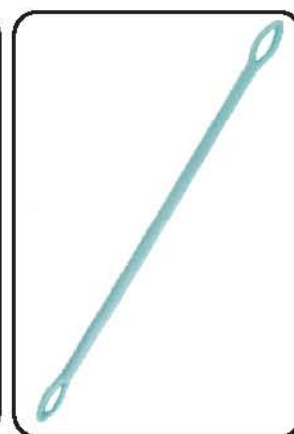

**Y2530 • Y2531**