

ORAL DOSING/REHYDRATION

DOSING/CROP NEEDLES


Stainless steel needles with an olive tip for easy dosing of birds and small mammals.

Straight needles. Length = 13cm.

Sizes:

- Y2521 12G
- /14 14G
- /16 16G
- /18 18G
- /20 20G

Shorter 8cm needles also available on request.

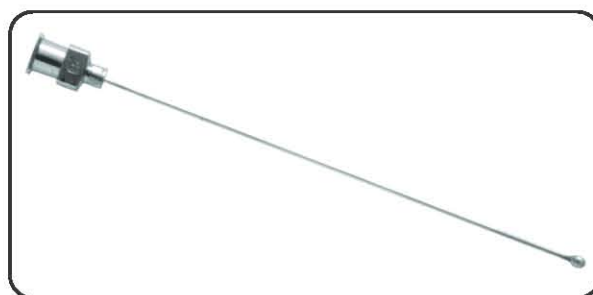

Y2521

DOSING/FLUSHING SYRINGES

- Y4415 A strong, 300ml, nylon/metal syringe with two different nozzles: the tapered (flushing) nozzle takes a uterine infusion pipette (Y3520) for easy flushing and the dosing nozzle enables easy dosing of oral medicaments.



Y4415

- Y4416 A 400ml nylon dosing syringe with the dosing nozzle to flush out horse's mouths during dental procedures



Y4416

PET NURSER KITS

Ideal for supplementing milk to newborn puppies and kittens (use with Milko-Pup or KittyMilk — each sold separately). Measurements are clearly marked on the side of the bottle and the kit includes a cleaning brush and extra silicone teats

Sizes:

- Y2620 50ml bottle
- Y2621 150ml bottle
- Y2620/1 Spare 50ml teats (5's)



Y2620 • Y2621

Y2535 PET PILL DISPENSER
To administer pills to small dogs and cats.



Y2535

Y2536 TABLET INTRODUCER (BUSTER)(10's)
A combination of a normal disposable syringe and an ordinary pet pill dispenser. The finger grips ensure full control during administration. Uses a small amount of water, and so utilises the normal swallowing reflex of the animal.



Y2536

STOMACH AND FEEDING TUBES

GASTRIC FEEDING TUBES

These sterile-packed silicone gastric feeding tubes have a bullet point with cap for luer fitting, graduated marking and four lateral eyes. Can be used long term.

Sizes:

- Y2590 6FR (2mm od x 40cm)
- Y2591 8FR (2,7mm od x 80cm)
- Y2592 12FR (4mm od x 120cm)
- Y2593 16FR (5,3mm od x 120cm)
- Y2594 20FR (7mm od x 120cm)



Y2590-Y2594

STOMACH TUBES, LARGE ANIMAL

Manufactured from non-toxic PVC tubing giving extra smooth finish, complete absence of tissue irritation. Resistance to kinking has been carefully considered in the tube design while retaining the necessary softness and flexibility. The tubes are fitted with a funnel end.

Sizes:

- Y4400 Foal (9,5mm od x 213cm)
- Y4401 Small (13mm od x 274cm)
- Y4402 Medium (16mm od x 274cm)
- Y4403 Large (19mm od x 304cm)



Y4400-Y4403

SECTION F – DOSING & REHYDRATION AIDS

STOMACH TUBES, K-9

PVC stomach tube with a lateral eye for bigger dogs. Useful in gastric lavage or administering larger volumes of fluids with a funnel.


Sizes:

Y4404 Medium (id 6.35mm) for dogs weighing greater than 18kg

Length: 110cm

Y4405 Large (id 7.87mm) for dogs weighing greater than 32kg

Length: 145cm

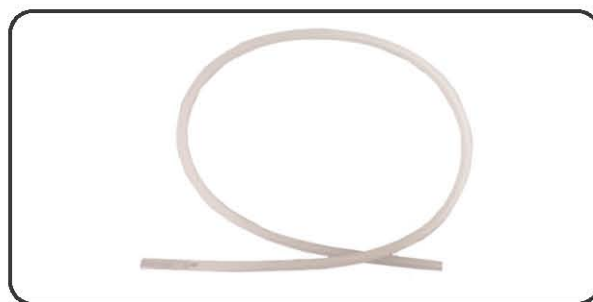


Y4404 • Y4405

IV ADMINISTRATION

Y4410 FLUTTER VALVE


Fast administration of IV fluids for bovines.

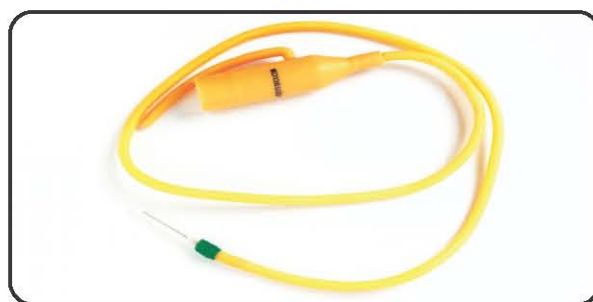

Y4410

Y4411 IV ADMINISTRATION SET, LARGE BORE

Increases flow rate up to 30L/hour. This volume is necessary in severe colic and critical care emergencies. The diameter is twice that of conventional IV drips; this translates into 16 times the flow rate of conventional sets.

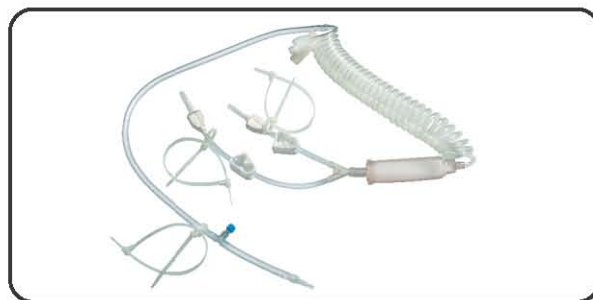


Y4411

Y4412 IV ADMINISTRATION EXTENSION TUBE, LARGE BORE

Size: 75cm

HI FLOW EQUINE IV CATHETER LARGE BORE (LONG TERM)

A large-bore catheter made of polyurethane, which softens at body temperature and can safely be left indwelling for up to 30 days. Some common catheter-related problems, such as phlebitis, can be minimised with these catheters. Non-thrombogenic and won't crack. Suture wings with holes.

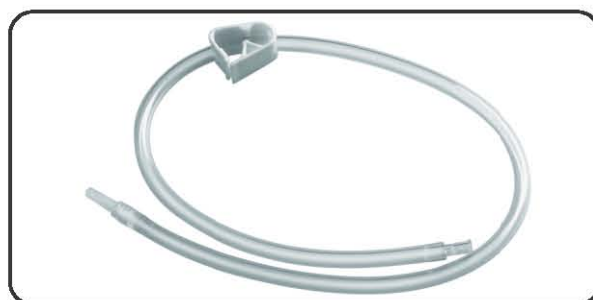

Y4412

SIZES (STERILE PACKS)

Y4413 14G (length = 8.75cm)

Y4414 12G (length = 13cm)

Y1098/I KYRON DUAL MODE INFUSION PUMP

Kyron's Dual Mode Infusion Pump is an intravenous infusion pump for various fluids, blood plasma extenders, chemotherapeutic agents, liquid nutrients, etc. It is applicable for micro and rapid infusion. The infusion pump has a user-friendly operation interface and bright LCD display which comes with 8 alarm levels for convenient operation in different environments.

The Infusion Pump comes standard with:

- Easy switch between drops/minute, ml/h and volume/time measurements
- Flow Rate Range of 1.0~1200 ml/h
- Bolus mode available
- Up to 20 IV sets can be added, stored and deleted in memory, convenient for operation after calibration
- Two year warranty

IV SPLINTS (SMALL ANIMALS)

COCHRAN IV SPLINT

While administering a drip, the cephalic vein is easily occluded when an animal is in sternal position with the elbow in flexion. The Cochran IV Splint allows the animal to lie comfortably yet still ensure a continuous IV flow.

Sizes:

Y3600 Small

Y3601 Medium

Y3602 Large

Y3603 Set of splints: small, medium, large & X-large



Y4413 • Y4414



Y1098/I



Y3600-Y3603