

CASTRATOR

BOVINE CASTRATOR, BLOODLESS
Castrators (Burdizzo-type) with cord-stop for bloodless castration.

Sizes:

- Y1130 23cm [can also be used for sheep]
- Y1131 30cm
- Y1132 40cm
- Y1133 47cm



Y1130-Y1133

- Y1103 **SHEEP CASTRATOR, ELASTRATOR**
Rubber ring applicator.



Y1103

EQUINE EMASCULATOR

- Y1100 **SERRA® EMASCULATOR**
A concave cutting plate draws the cord concentrically into the jaws of the castrator, so avoiding slipping of the cord. The cord and surrounding strands are centrally gripped, pressed together and cut through. Reduces tissue slippage and post-op haemorrhage.


Length: 35cm



Y1100

- Y1101 **RIEMERS EMASCULATOR**
This emasculator has a double crushing part with a separate cutting lever, which is worked separately. This allows extra time for crushing action and eliminates after-bleeding.

Length: 31 cm

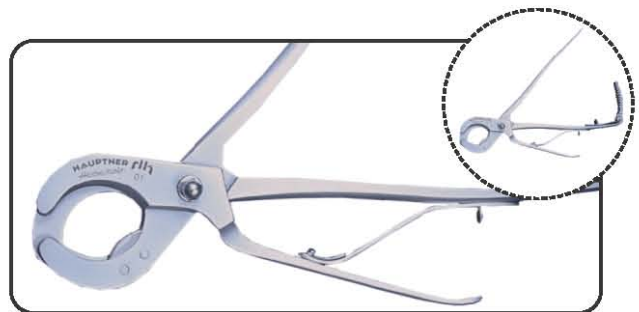

Y1101

- Y1102 **SAND EMASCULATOR**
Twin- and single-blade for crushing only.



Y1102