

**ARTERY/HAEMOSTATS**

*Self-locking stainless steel artery/haemostatic forceps are used to grasp the artery/blood vessels during ligation in order to prevent bleeding.*


**FERGUSON ANGIOTRIBE FORCEPS**

Heavy-duty forceps with a male/female groove down the length. Ideal for grasping fallopian tubes of large dogs during spays; as well as spermatic cords during castrations.

*Length: 19cm*

Y1701 Straight

Y1700 Curved

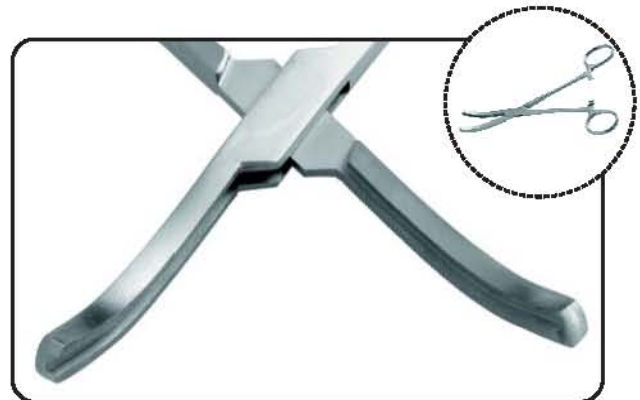

Y1700


**HALSTEAD MOSQUITO ARTERY FORCEPS**

Precision, fine-pointed forceps, especially useful for handling small vessels.

*Length: 12cm*

Y1710 Straight

Y1711 Curved

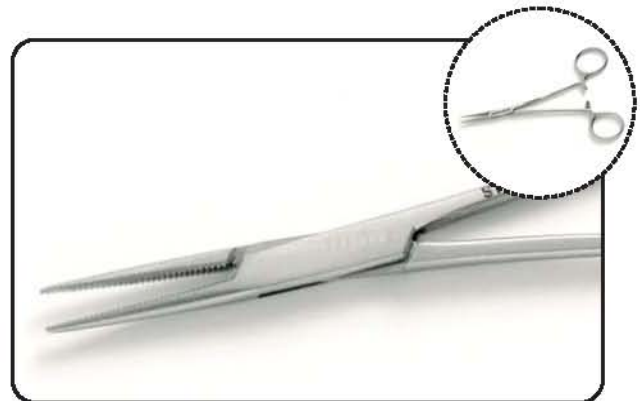

Y1710

**SPENCER WELLS ARTERY FORCEPS**

Straight:

Size:

Y1720 13cm

Y1722 15cm

Y1724 18cm


Curved:

Size:


Y1721 13cm

Y1723 15cm

Y1725 18cm

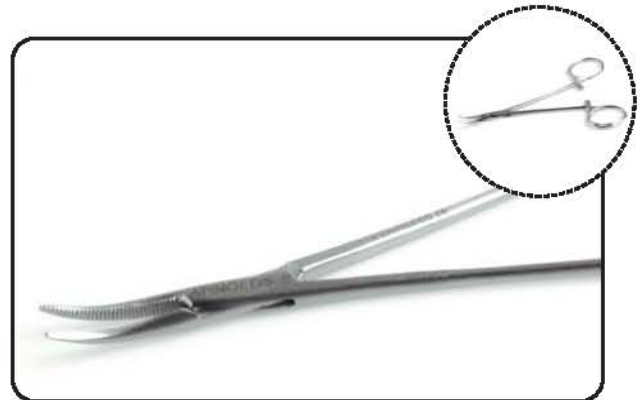

Y1711

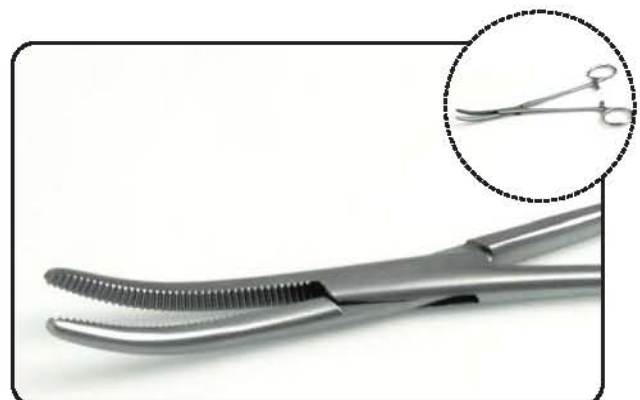

Y1721 • Y1723 • Y1725

## DISSECTING

Non-locking stainless steel “thumb” forceps are used to grasp delicate tissue during fine dissection work. Available either as serrated (with small serrations at the tip for dressings, etc) or rattooth (with a single “tooth” on one side of the tip, which fits into a groove on the opposing side of the tip and useful for holding tissue).

### Y1500 ADSON BROWN DISSECTING FORCEPS

Have fine micro-teeth (8/9 teeth) down the length of the tip for a more atraumatic grip (e.g. skin).  
Length: 12cm

ADSON DISSECTING FORCEPS  
Length: 12cm

Y1501 Rattooth

Y1502 Serrated

### Y1509 JEAN'S DISSECTING FORCEPS

Heavy-duty forcep with a rattooth tip and some serrations on the shaft directly behind the tip (ideal for cattle work).  
Length: 14cm Rattooth

### STANDARD DISSECTING FORCEPS

Serrated (plain/dressing forceps):

Sizes:

Y1503 13cm

Y1505 15cm

Y1507 18cm


Rattooth (tissue forceps):

Sizes:


Y1504 13cm

Y1506 15cm


Y1508 18cm

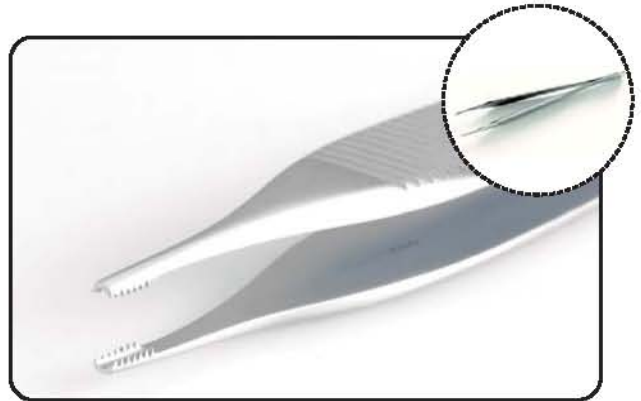


Y1500

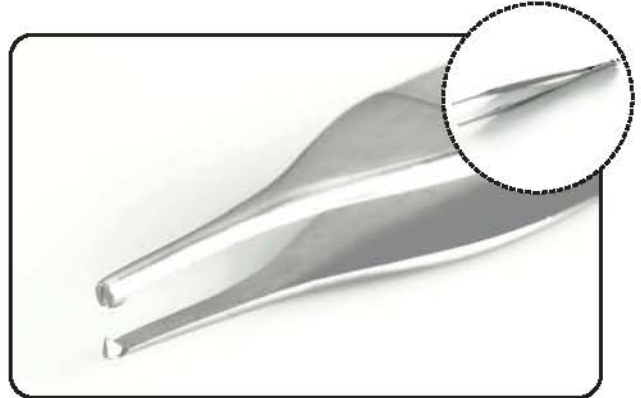

Y1501

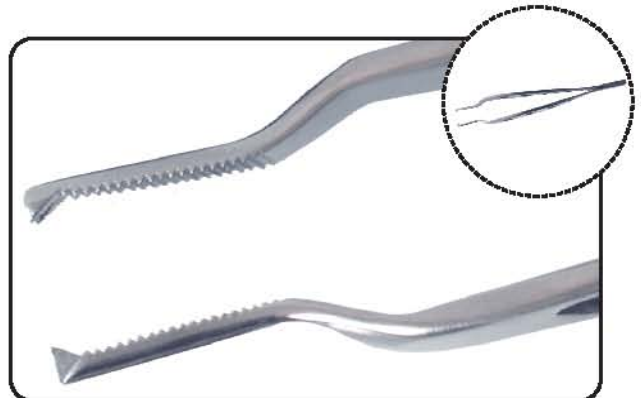

Y1509

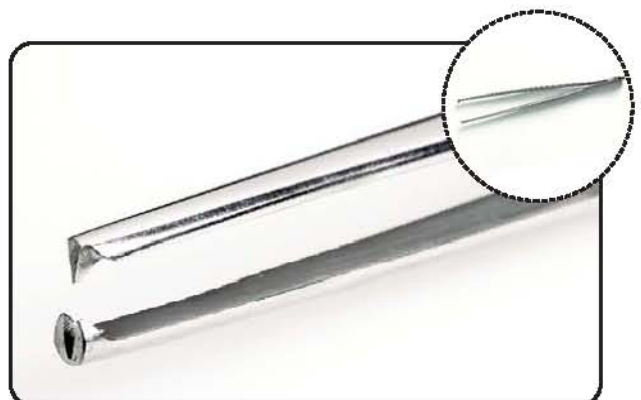

Y1504 • Y1506 • Y1508

**EAR**

Non-locking aural forceps are stainless steel with a lengthy shaft and a hinged tip. Useful for removal of foreign bodies, wax, etc from an enclosed space (i.e. ear canal).

**HARTMAN AURAL FORCEPS**

Sizes (shaft length):

Y1650 9cm

Y1651 13cm

Y1652 16,5cm

20cm available on request.

**Y1660 JOBSON HORNE'S PROBE**

A very useful instrument, technically not a forceps, used for removal of wax and cleaning of ears. One end functions as an applicator: with serrated edges it is designed to hold cotton-wool. The other end is a blunt curette for removal of ear wax and cleaning of ears.

Length: 18cm

**EYE**

"Thumb-style" stainless steel forceps for fine dissection during ocular surgery.

**IRIS DISSECTING FORCEPS**


Length: 11cm

Y1511 Rattooth


Y1510 Serrated

**Y1617 INSTRUMENT PIN (MAYO)**


An easy way to keep instruments together in a pack.

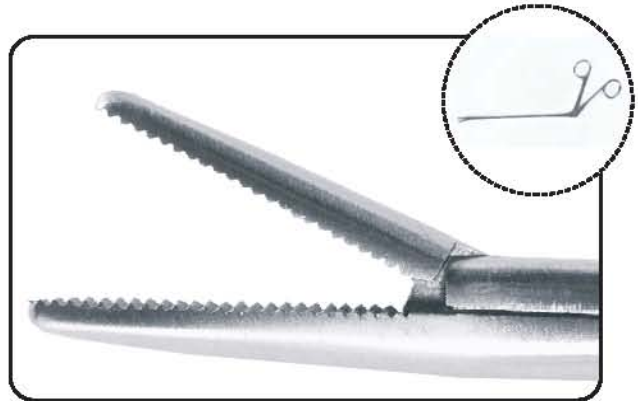


Y1652

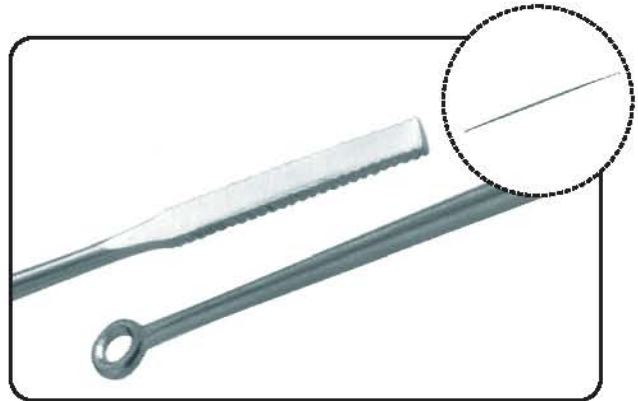

Y1660

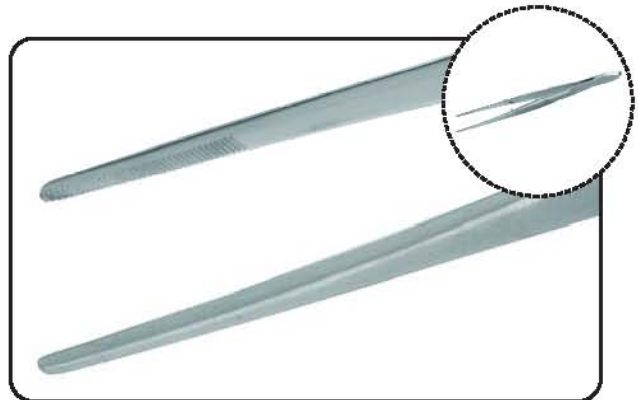


Y1510

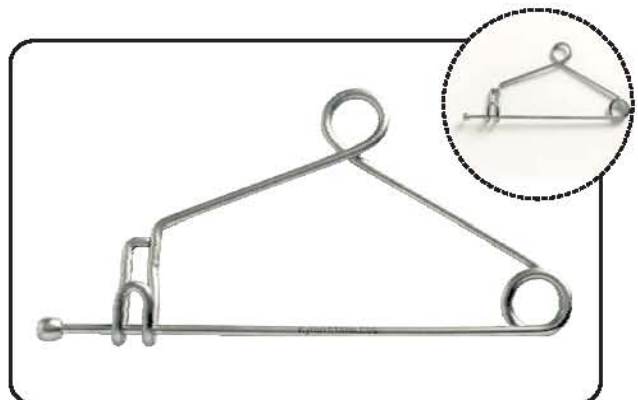

Y1617

## INTESTINAL

Self-locking stainless steel intestinal forceps, which are non-crushing and which contain serrations down the length of the shaft. Used to temporarily occlude sections of bowel.

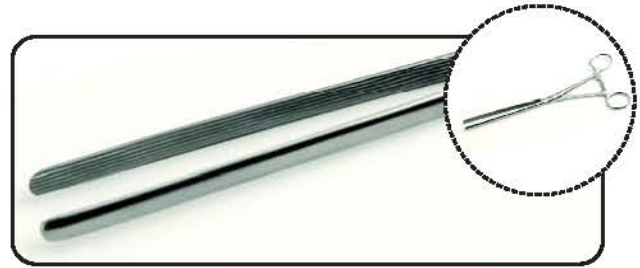


Y1780

## DOYENS INTESTINAL FORCEPS

Length: 23cm

- Y1780 Straight
- Y1781 Curved

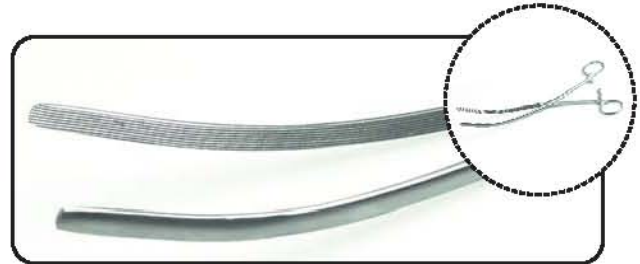


Y1781

## Y6031 OESOPHAGEAL FOREIGN BODY FORCEPS

An especially long forceps that can remove foreign bodies such as fish hooks or lodged bones from even the distal oesophagus or trachea. The shaft diameter of 4.2mm is as small as possible to allow the placement of an endoscope alongside. In addition, the backward facing teeth allow for maximum grip, which reduces the chance of slipping.

Shaft Length: 50cm (4.2mm diameter)

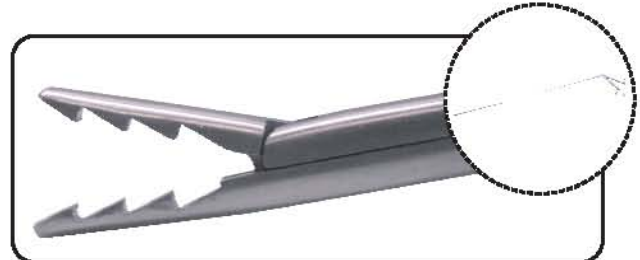

Y6031

## TISSUE

Self-locking stainless steel tissue forceps with interdigitating teeth to grip tissue.

- Y1600 ALLIS TISSUE FORCEPS (4/5 teeth)  
Length: 16cm  
(19cm available on request)



Y1600

## TOWEL CLIPS


### BACKHAUS TOWEL CLIPS

Sizes:

- Y1615 9cm
- Y1616 13cm

### Y1610 X-ACTION TOWEL CLIPS

Size: 9cm

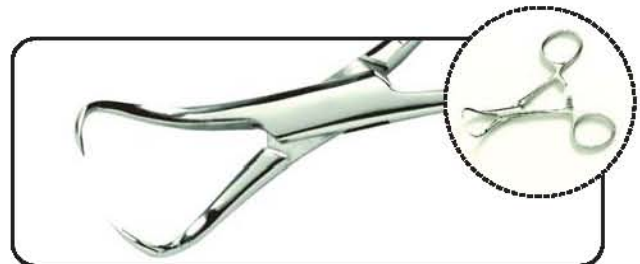


Y1615



Y1610

**WOUND RETRACTORS:  
HAND-HELD**

Y1635 DESMARRE'S EYELID RETRACTOR  
*Length: 13,5cm - Width: 1cm*

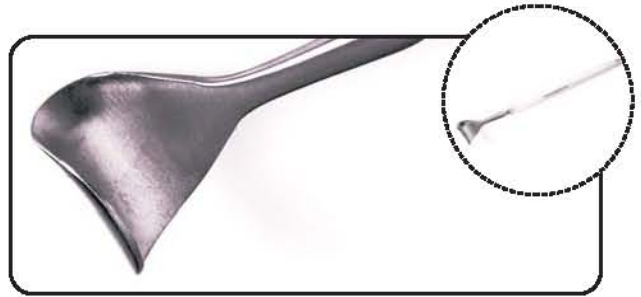


Y1635

Y1627 SENN MILLER RETRACTOR  
*Size: 17cm*



Y1627


Y1630 SPAY HOOK (blunt)  
*Length: 16cm*

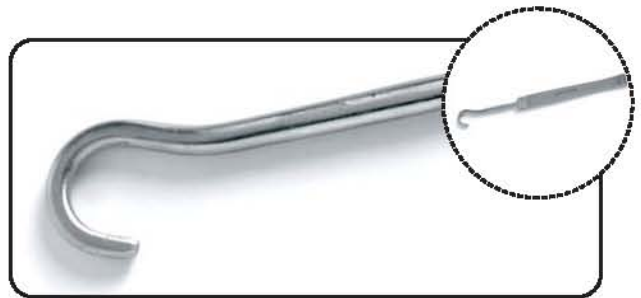

Y1630

**WOUND RETRACTORS:  
SELF-RETAINING**

Y1626 GELPI SELF-RETAINING  
RETRACTOR  
*Size: 18cm*

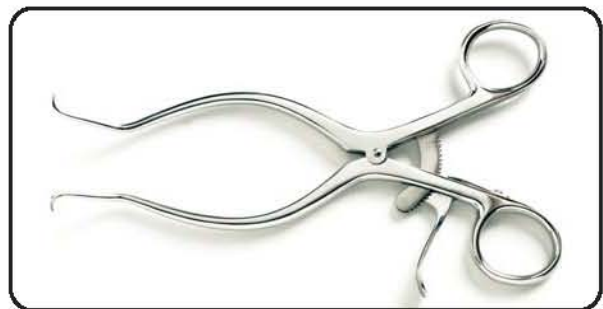

Y1626

WEITLANER SELF-RETAINING  
RETRACTOR

*Sizes:*

Y1624 16cm

Y1625 20cm



Y1624 • Y1625