


SCALPEL HANDLES & BLADES


BARBER OF SHEFFIELD

Barber of Sheffield scalpel blades are made of the finest quality steel to the most stringent manufacturing quality control standards. The cutting edge is of the highest possible standard. The blades are hermetically sealed in individual aluminium foil easy-to-peel sachets and are sterilised by gamma radiation. 100 Sterile blades per box.

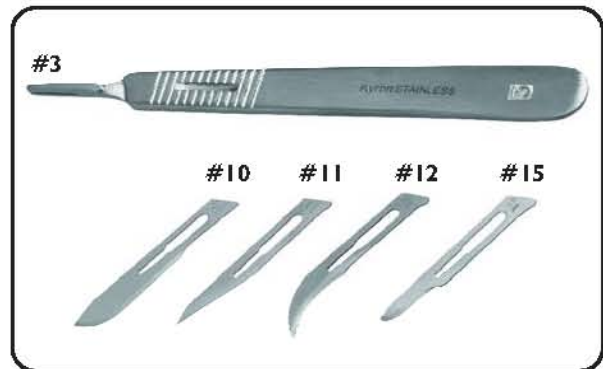

Y3060 • Y4509-Y4512

- Y3060 Handle # 3 to fit blades:
- Y4510 #10
- Y4509 #11
- Y4511 #12
- Y4512 #15

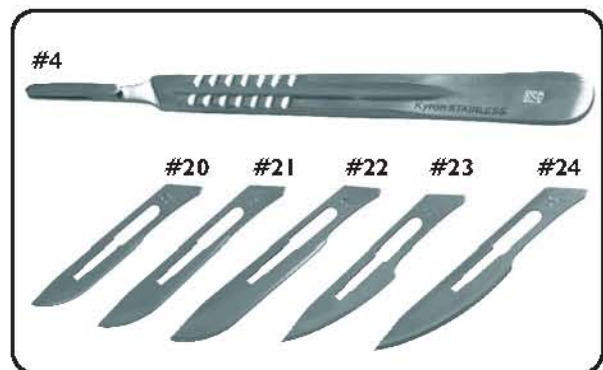

Y3061 • Y4513-Y4517

- Y3061 Handle # 4 to fit blades:
- Y4513 #20
- Y4514 #21
- Y4515 #22
- Y4516 #23
- Y4517 #24

POST MORTEM

- Y3200 Handle # 40 to fit blade:
- Y3201 # 40 (10's)



Y3200 • Y3201


SCISSORS

IRIS SCISSORS

Initially named so because of their usefulness in eye surgery but now widely used, Iris Scissors are small, delicate (but very sharp) scissors with a pointed end — often used for fine dissection work.

Sizes:

- Y3400 Straight (11cm)
- Y3401 Curved (11cm)

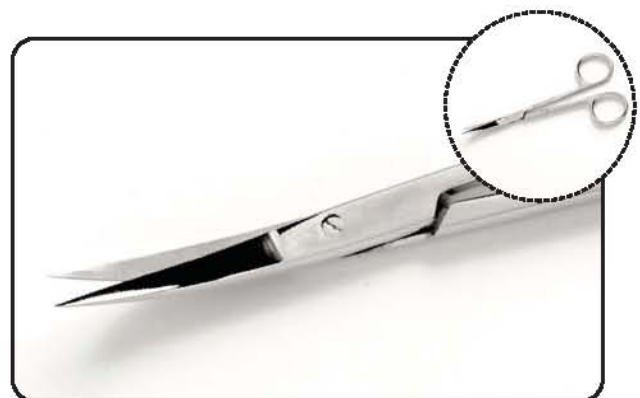

Y3401

MAYO SCISSORS

A more heavy-duty type of scissors with semi-blunt, pointed blades. Mayo Scissors are heavier than Metzenbaum Scissors and can be used to cut fascia and through tissue.

STRAIGHT (ECONOMY):

Sizes:

Y3425 14cm

Y3427 17cm

STRAIGHT (GIMMI — German quality scissors with sharp, precisely engineered cutting edges):

Sizes:

Y3430 14cm

Y3432 17cm

CURVED (ECONOMY):

Sizes:

Y3426 14cm

Y3428 17cm

CURVED (GIMMI):

Sizes:

Y3431 14cm

Y3433 17cm

METZENBAUM SCISSORS

Have a narrow, rounded tip and are more delicate than Mayo Scissors (and have a longer handle:blade ratio than the Mayo). Metzenbaum Scissors are used to cut or dissect delicate tissue (avoid using these scissors to cut sutures as this will dull the blades).

STRAIGHT (GIMMI):

Sizes:

Y3460 14,5cm


Y3462 18cm

CURVED (GIMMI):

Sizes:

Y3461 14,5cm


Y3463 18cm

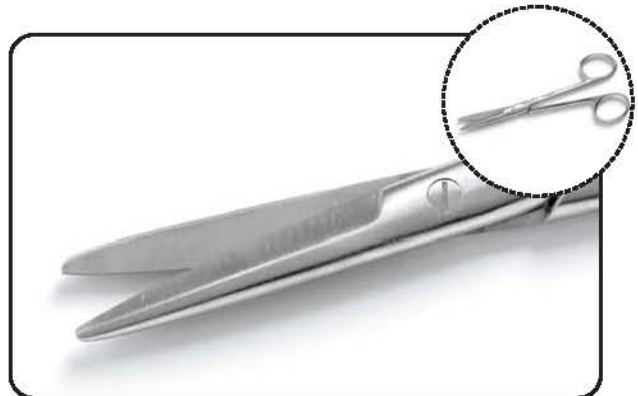


Y3425 • Y3427 • Y3430 • Y3432



Y3426 • Y3428 • Y3431 • Y3433

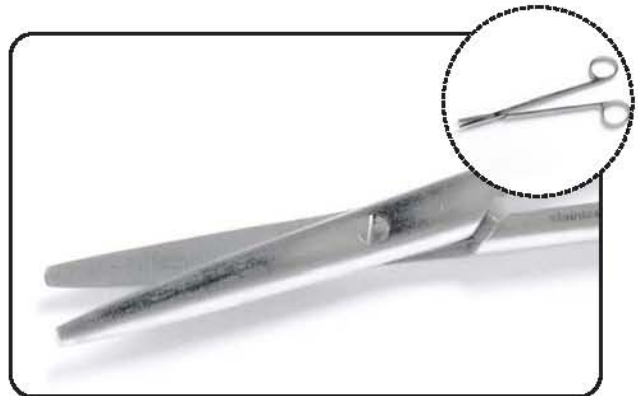

Y3460 • Y3462

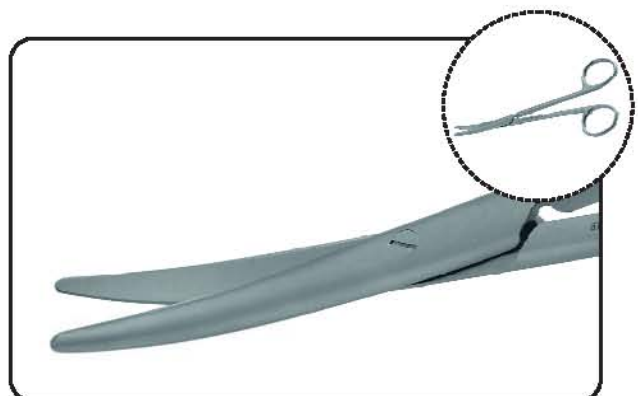

Y3461 • Y3463

STANDARD OPERATING SCISSORS

May be used for e.g. cutting suture material or cutting hair on the leg when preparing a vein for an IV.

SHARP/SHARP STRAIGHT:

Sizes:

Y3411 13cm

Y3414 15cm

SHARP/BLUNT STRAIGHT:

Sizes:

Y3410 13cm

Y3413 15cm

BLUNT/BLUNT STRAIGHT:

Sizes:

Y3409 13cm

Y3415 15cm

BLUNT/BLUNT CURVED:

Sizes:

Y3412 13cm

Y3416 15cm

STITCH-REMOVAL SCISSORS

Scissors with a modified tip for easy removal of stitches by inserting the one tip beneath the suture.

Y3437 HEATH'S STITCH-REMOVAL SCISSORS

Size: 16cm

SPENCER STITCH-REMOVAL SCISSORS

Sizes:

Y3435 9cm

Y3436 11cm

Y3440 UNIVERSAL SCISSORS

A very useful multi-purpose scissor (handy in the home, too). Size: 19cm

Y3445 VESSEL SCISSORS

Ideal for urethrostomies — blunt tip prevents accidental penetration of the urethra. Size: 11,5cm



Y3411 • Y3414



Y3410 • Y3413



Y3412 • Y3416



Y3437



Y3435 • Y3436



Y3440



Y3445