

NEEDLES

ANCHOR SUTURE NEEDLES

Made in the USA, finest quality re-usable stainless steel suture needles with reverse cutting edge. Reverse cutting needles have the cutting edge on the outside and the flat surface towards the wound and this may offer better resistance against tissue tear. Reverse cutting needles however do require more force to penetrate the tissue. Packed in 6's. Actual sizes illustrated.

$\frac{1}{2}$ CURVED, CUTTING


- Y2710 # 1
- Y2711 # 2
- Y2712 # 4
- Y2713 # 6
- Y2714 # 8
- Y2715 # 10
- Y2716 # 12
- Y2717 # 14
- Y2718 # 16
- Y2719 # 18
- Y2720 # 20

CURVED ($\frac{3}{4}$ CIRCLE), CUTTING


- Y2730 # 1
- Y2731 # 2
- Y2732 # 4
- Y2733 # 6
- Y2734 # 8
- Y2735 # 10
- Y2736 # 12
- Y2737 # 14
- Y2738 # 16
- Y2739 # 18
- Y2740 # 20

Y2746 PM CUTTING


Double curved needle: 12cm.
Illustration is not actual size.

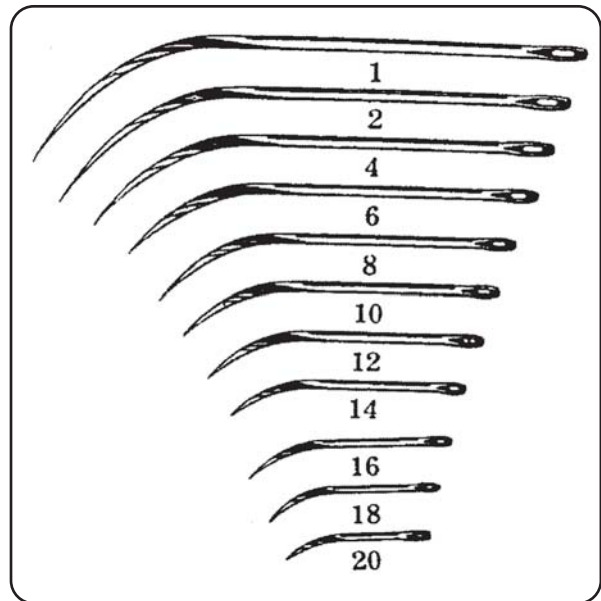

Y2710-Y2720 (actual sizes)

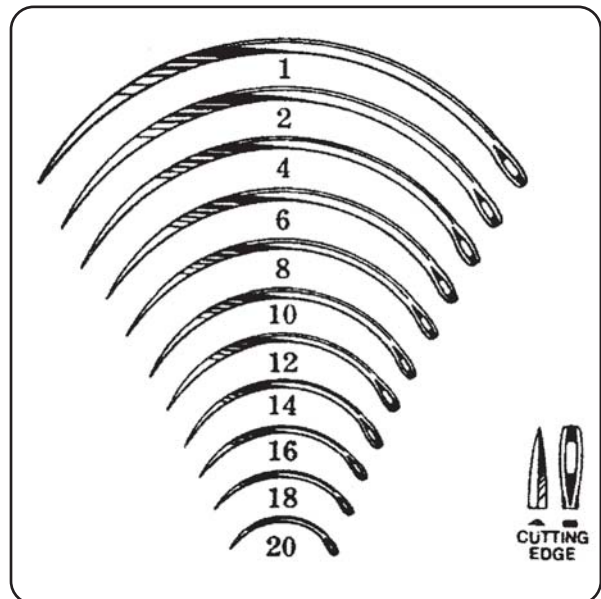

Y2730-Y2740 (actual sizes)

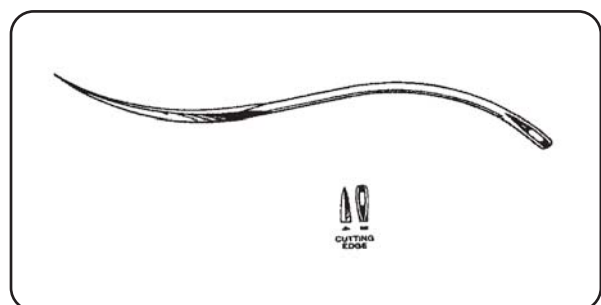

Y2746 (not actual size)

BARBER of SHEFFIELD SUTURE NEEDLES

Finest quality re-usable stainless steel suture needles with triangular cutting edge or round/taper point.

Conventional cutting needles have the cutting edge on the inside of the needle curvature and hence directed towards the wound. Suture tension is therefore at the top of the triangular part of the needle and tissue tear resistance is weak. Made in the UK. Packed in 12's. *Actual sizes illustrated.*

½ CURVED, CUTTING

Y5020	# 1
Y5021	# 2
Y5022	# 4
Y5023	# 6
Y5024	# 8
Y5025	# 10
Y5026	# 12
Y5027	# 14
Y5028	# 16
Y5029	# 18
Y5030	# 20




CURVED TRIANGULAR ($\frac{3}{8}$ CIRCLE), CUTTING

Y5000	# 1
Y5001	# 2
Y5002	# 4
Y5003	# 6
Y5004	# 8
Y5005	# 10
Y5006	# 12
Y5007	# 14
Y5008	# 16

CURVED ($\frac{3}{8}$ CIRCLE), ROUND BODIED / TAPER

Round bodied needles taper gradually to a point whereas triangular bodies have cutting edges along three sides.

Y5040	# 1
Y5041	# 2
Y5042	# 4
Y5043	# 6
Y5044	# 8
Y5045	# 10
Y5046	# 12
Y5047	# 14
Y5048	# 16
Y5049	# 18
Y5050	# 20

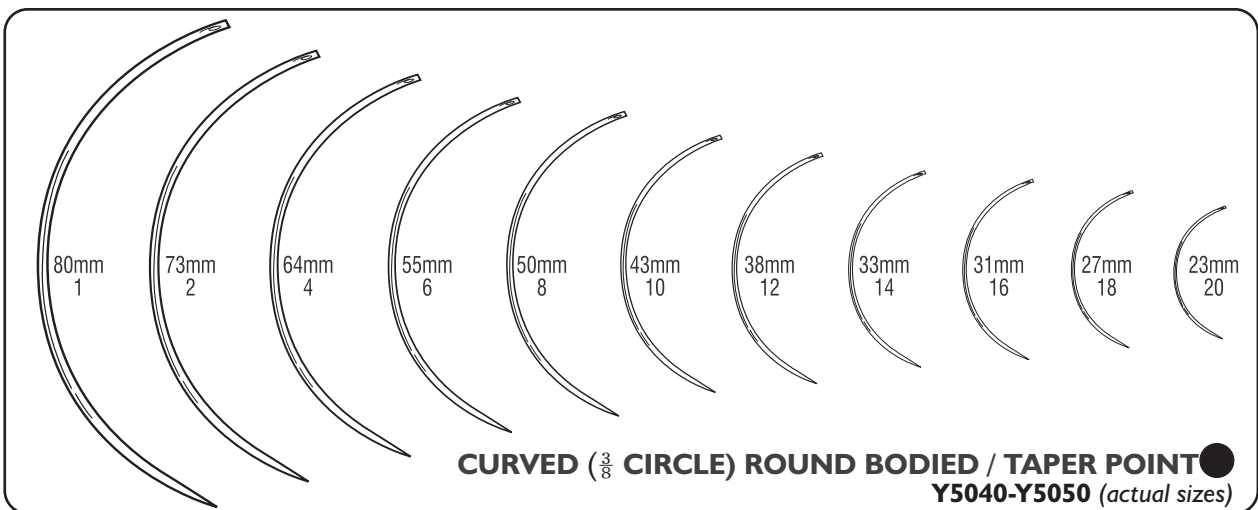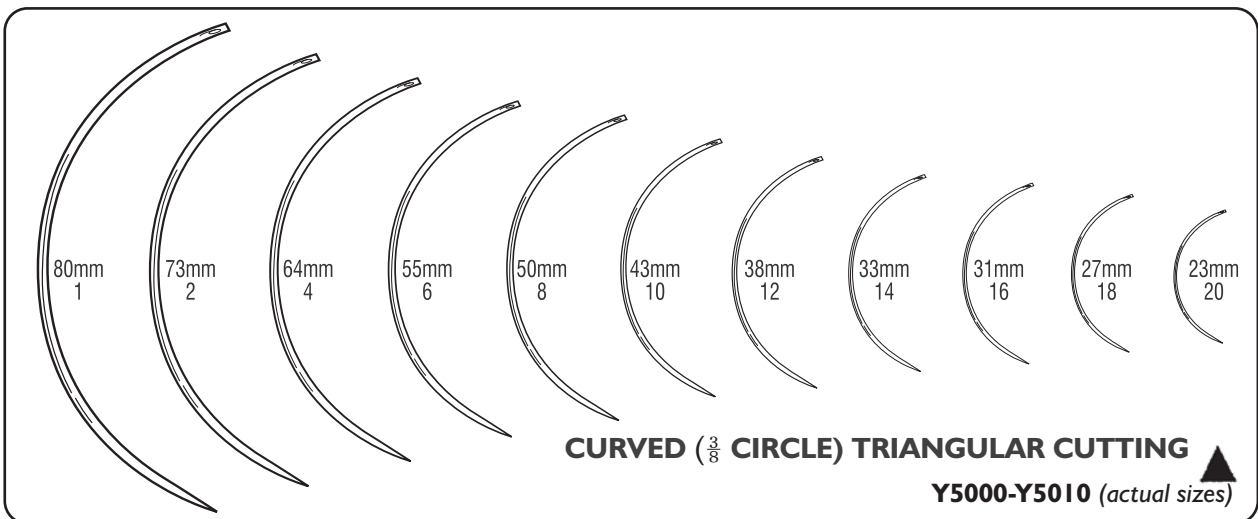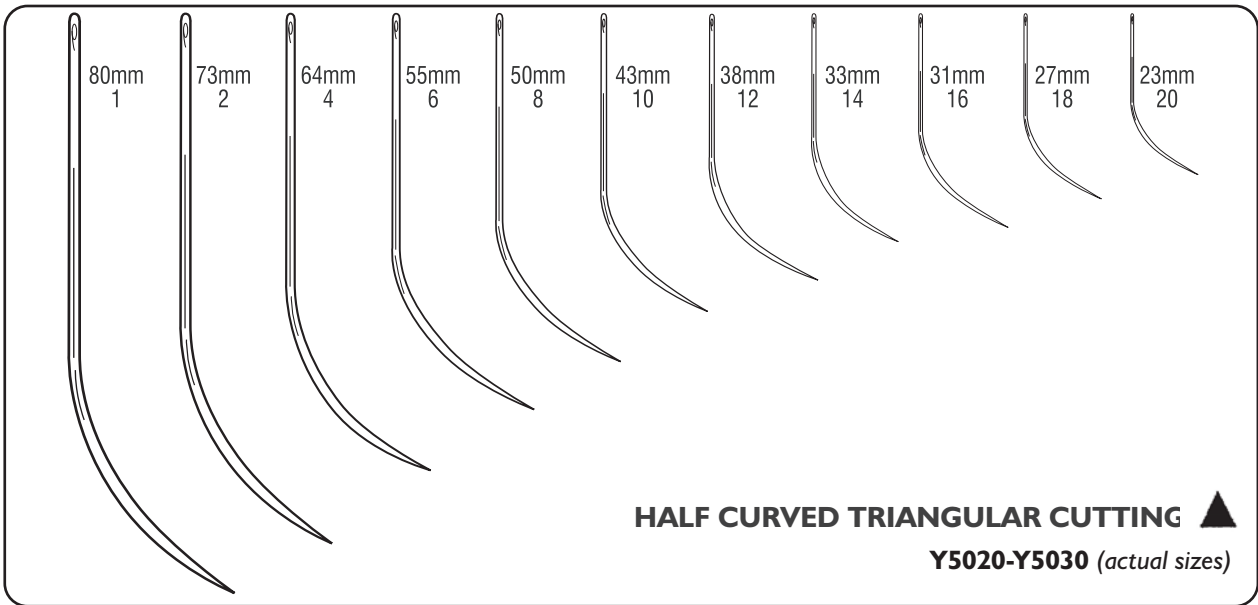


SECTION V — SUTURING

STRAIGHT CUTTING

- Y2790 # 2 (7cm)
- Y2791 # 4 (6cm)
- Y2792 # 8 (5cm)
- Y2793 # 14 (4cm)

Sharps container for the disposal of used suture and surgical needles.

- Y2794 4L container
- Y2795 0.5 L container

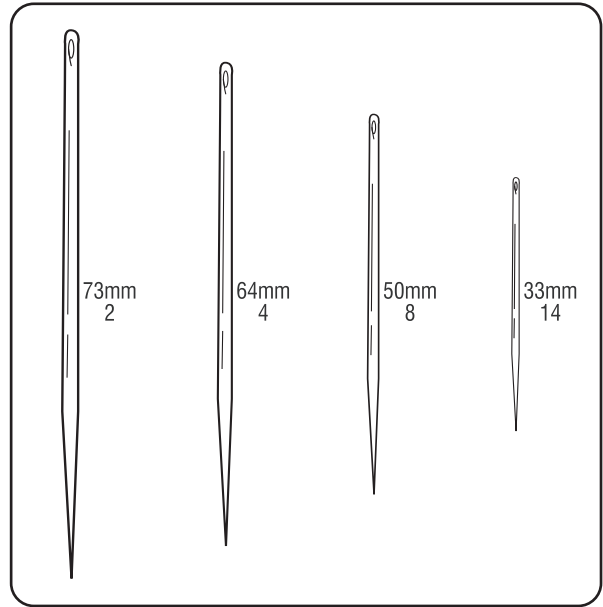

Y2790-Y2793 (actual sizes)



Y2794-Y2795