



Y4202 TEAT BLADES

These stainless steel blades are designed with a specially shaped cutting edge inside a hook, making a cut in the sphincter so that the outside opening of the teat canal remains whole. Packed in 6's.

Length: 5cm

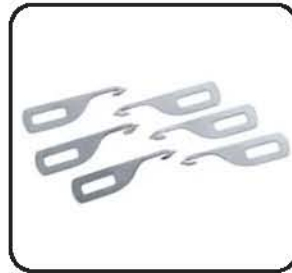

Y4202

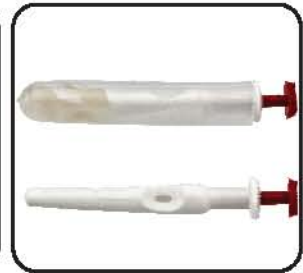

Y4225

TEAT CANNULA

Y4225 BOVIVET TEAT DILATOR PLUGS (10's)



Y4221

Y4221 UDDER INFUSION CANNULA

This stainless steel cannula with a luer adapter is used for flushing teats or infusing medication into teats. Also used for abdominal taps in horses.

Length: 78mm



Y4205

Y4205 TEAT DILATOR

Length: 11,8cm



Y4206

Y4206 TEAT DILATOR WITH SCREW

Length: 16,8cm

TEAT KNIFE

Y4203 FRENCH PATTERN TEAT KNIFE

Length: 10,5cm



Y4203

Y4200 HUG'S TEAT KNIFE/LANCET

Length: 15cm

Blade: 15mm




Y4200

SECTION W — TEAT INSTRUMENTS

Y4201 TEAT SLITTER, DANISH PATTERN

Blades are opened by loosening the small burr on the side of the instrument.

Length: 19cm

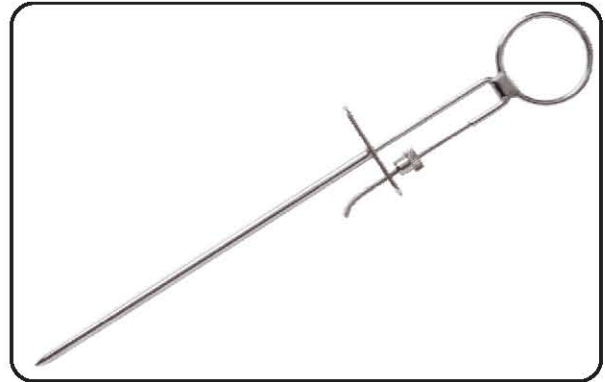

Y4201

Y4222 TEAT SYPHONS, HAUPTNER

Cannula with bulbous end to take rubber tubing.


Size: 8,5cm



Y4222

Y4210 TEAT TUMOUR EXTRACTOR

With two cutting edges that come together for removal of fibrous tissue.

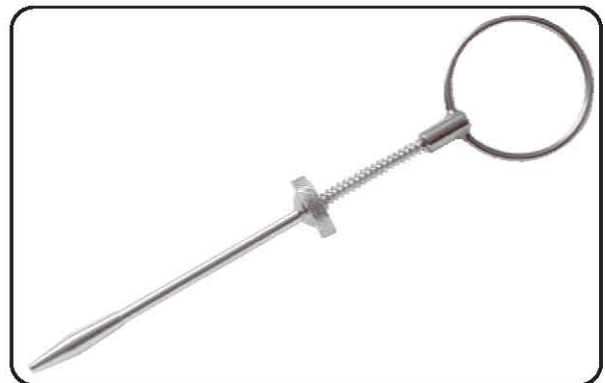

Y4210